

## Occult blood in feces 2024.1



<p><b>Scheme Period</b></p>	<p>Occult blood in feces 2024.1 January 1, 2024 - March 31, 2024</p>
<p><b>Report for</b></p>	
<p><b>Subscriptions Supervision and report authorisation</b></p>	<p>76 dr. H. van der Vuurst (Coördinator) reachable at office@skml.nl</p>
<p><b>Reference</b></p>	<p><a href="#">MUSE manual</a> and <a href="#">Analyte specifications</a></p> <p>On the SKML website you will find the MUSE manual and an overview of all analytes of all schemes. For each analyte the overview shows the source of the target value, the tolerance and what it is based on. The type of material and the specification of the commutability, homogeneity and stability have been included as well. Older versions of these documents are available upon request.</p>

<p><b>Qualitative</b></p>	<p><b>reported</b></p> <p>100%</p>
---------------------------	------------------------------------



### Comments

Monsters zijn gespiked met: 80 (2024.1A), 99,6 (2024.1B), 0 (2024.1C), 300 (2024.1D), 1200 (2024.1E), 0 (2024.1F) ug Hb/L.

## Occult blood in feces 2024.1

Analyte	Unit	Trueness				Precision		Performance			
		your mean	ref.	cons.	SDbl	your SD	SDwl	this survey	PS	cumulative	PSc
Hb in feces	ug/L	341		307	79	17	31				


### Qualitative scores



Analyte	This survey				Cumulative			
	correct	incorrect	total	pictogram	correct	incorrect	total	pictogram
Hb in feces (qual.)	6	0	6		6	0	6	


# Occult blood in feces 2024.1

Hb in feces

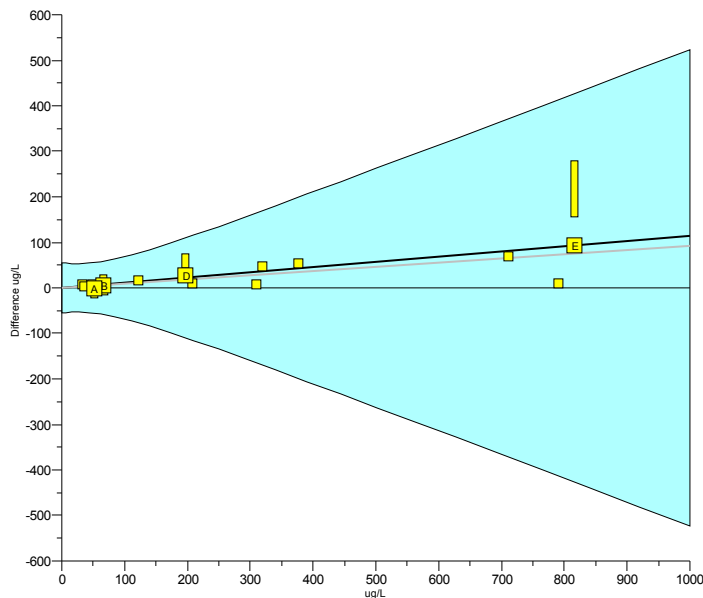units : ug/L



	2024.1	cumulative
Trueness	+11%	+9.1%
Precision	1.6%	6.2%
Number	4	16
Outliers	0	0
Sigma-TE		
Sigma-SA	6.0	6.0
Score pictogram		
Regression line	$0 + 1.114.x$	$0 + 1.092.x$
Consensus group	Sentinel method W	
Method	Sentinel method W	
Analyser	Sentifit 270	
Your factor	$0.0 + 1.000.x$	
Method factor	$0.0 + 1.022.x$	

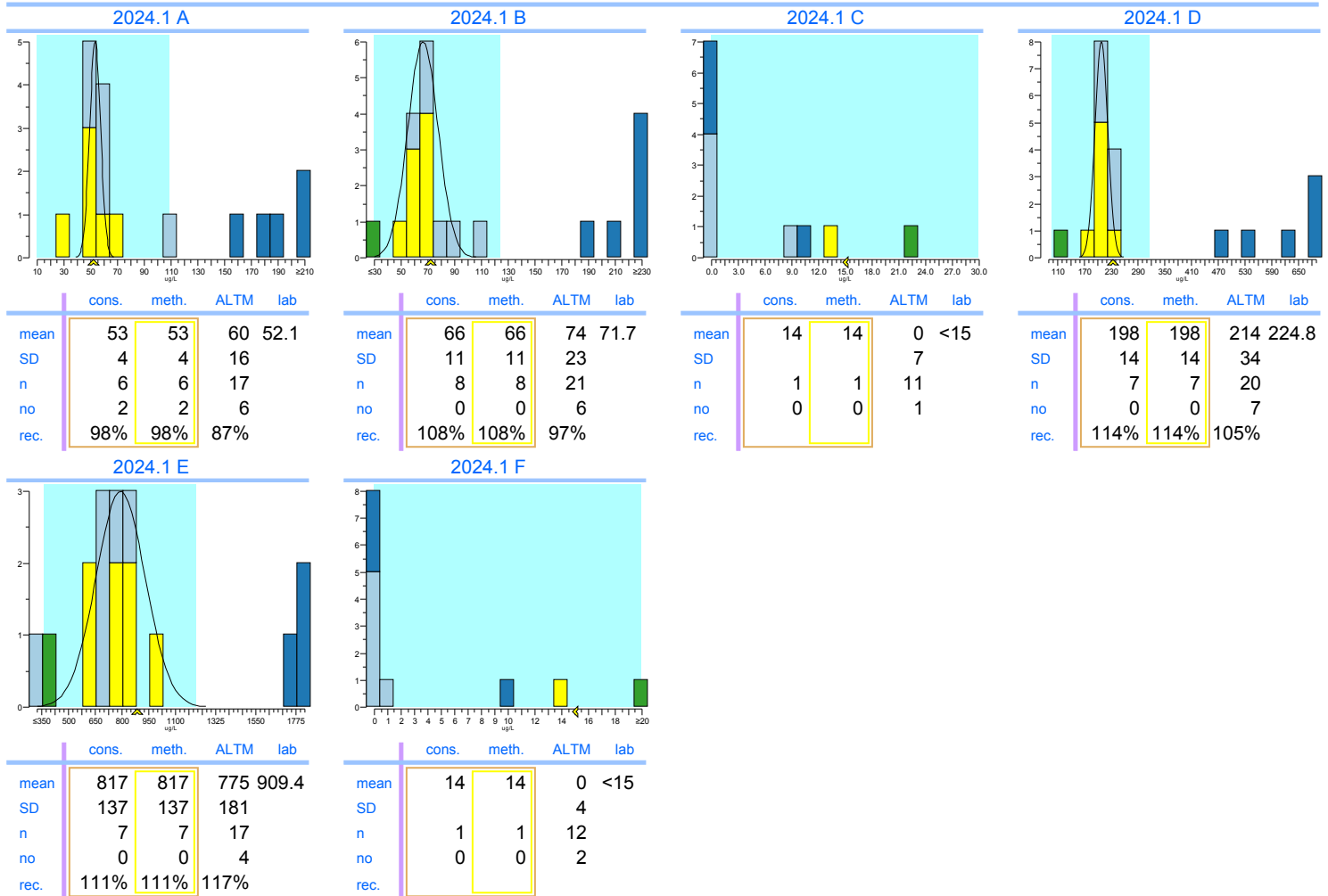

Legend




# Occult blood in feces 2024.1

Hb in feces (qual.)


units : qualit.

2024.1 A

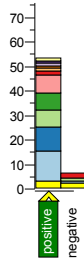


score : 2

2024.1 E

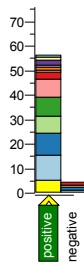


score : 2

2024.1 B

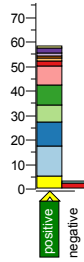


score : 2

2024.1 F

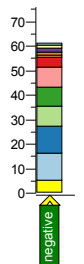


score : 2

2024.1 C

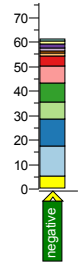

score : 2

2024.1 D

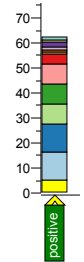

score : 2

Legend

- |                                                                                     |                                                                               |                                                                                        |                                                                           |                                                                                   |
|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------|----------------------------------------------------------------------------------------|---------------------------------------------------------------------------|-----------------------------------------------------------------------------------|
| <span style="border: 1px solid yellow; padding: 2px;">■ Sentinel</span>             | <span style="border: 1px solid lightblue; padding: 2px;">■ Alere Actim</span> | <span style="border: 1px solid blue; padding: 2px;">■ Innovacon</span>                 | <span style="border: 1px solid lightgreen; padding: 2px;">■ Others</span> | <span style="border: 1px solid green; padding: 2px;">■ Human</span>               |
| <span style="border: 1px solid pink; padding: 2px;">■ Medix</span>                  | <span style="border: 1px solid red; padding: 2px;">■ Eiken</span>             | <span style="border: 1px solid orange; padding: 2px;">■ DA/Cortez Diagnostics</span>   | <span style="border: 1px solid orange; padding: 2px;">■ Axon Lab</span>   | <span style="border: 1px solid lightgrey; padding: 2px;">■ Alere Clearview</span> |
| <span style="border: 1px solid purple; padding: 2px;">■ Standard Diagnostics</span> | <span style="border: 1px solid yellow; padding: 2px;">■ Manual</span>         | <span style="border: 1px solid teal; padding: 2px;">■ CerTest FOB one step card</span> |                                                                           |                                                                                   |

# Occult blood in feces 2024.1

## Time plots

Hb in feces

ug/L

External

